


BioVendor

Research
and Diagnostic Products

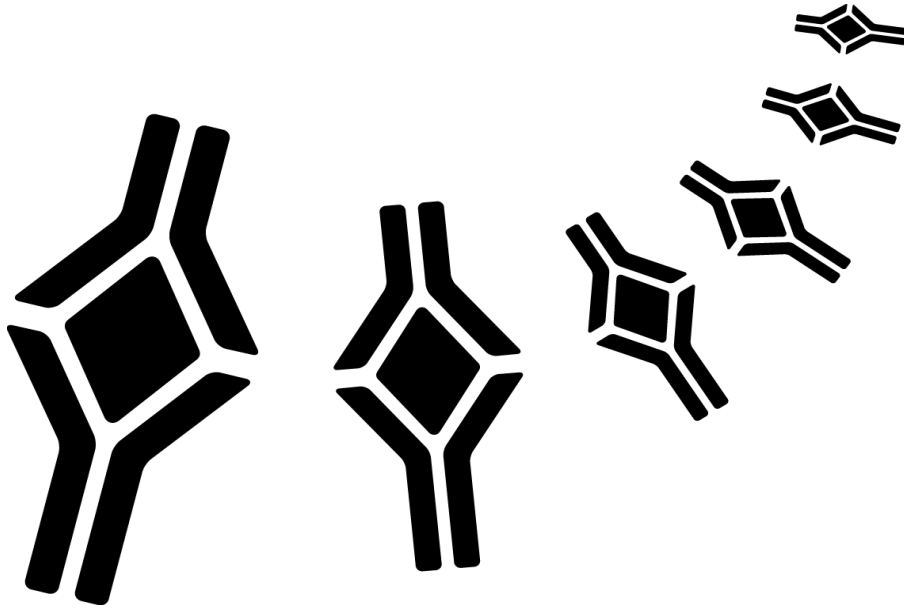

HUMAN CRP ELISA

Product Data Sheet

Cat No.: RAP001

European
Union:



Rest of the world:
For research use only!

CONTENTS

1.	INTENDED USE	3
2.	INTRODUCTION	3
3.	PRINCIPLE OF THE CRP ELISA	4
4.	REAGENTS	4
5.	MATERIALS REQUIRED BUT NOT SUPPLIED	5
6.	WARNINGS AND PRECAUTIONS FOR USERS	5
7.	STORAGE CONDITIONS	5
8.	SPECIMEN COLLECTION AND PREPARATION	5
9.	ASSAY PROCEDURE	6
10.	CALCULATIONS	8
11.	PERFORMANCE CHARACTERISTICS	9
12.	CLINICAL DATA	9
13.	TROUBLESHOOTING	10
14.	REFERENCES	10

**➤➤ This kit is manufactured by:
BioVendor – Laboratorní medicína a.s.**

➤➤ Use only the current version of Product Data Sheet enclosed with the kit!

1. INTENDED USE

Enzyme Immunoassay for the Quantitative Determination of CRP in Human Serum and Plasma.

2. INTRODUCTION

C-Reactive Protein (CRP) is an acute-phase protein, produced exclusively in the liver. Interleukin-6 is the mediator for the synthesis by the hepatocytes of CRP, a pentamer of approximately 120,000 Daltons. CRP is present in the serum of normal persons at concentrations ranging up to 5 mg/l. The protein is produced by the foetus and the neonate and it does not pass the placental barrier, as such it can be used for the early detection of neonatal sepsis.

Because febrile phenomena, leukocyte count and erythrocyte sedimentation rate (ESR) are often misleading, investigators and clinicians now prefer a quantitative CRP determination as a marker for acute inflammation and tissue necrosis. Within 6 hours of an acute inflammatory challenge the CRP level starts to rise.

Serum concentration of CRP increases significantly in cases of both infectious and non-infectious inflammation, of tissue damage and necrosis and in the presence of malignant tumours. CRP is present in the active stages of inflammatory disorders like rheumatoid arthritis, ankylosing spondylitis, Reiter's syndrome, psoriatic arthropathy, systemic lupus erythematosus, polyarteritis, ulcerative colitis and Crohn's disease.

Injuries causing tissue breakdown and necrosis are associated with increases in serum CRP which are seen in thermal burns, major surgery and myocardial infarction.

Widespread malignant disease with carcinoma of the lung, stomach, colon, breast, prostate and pancreas, Hodgkin's disease, non-Hodgkin's lymphoma and lymphosarcoma will give rise to high levels of CRP resulting from tissue damage by invading tumour cells. Therefore, CRP may be used to monitor malignancy.

The CRP-level increases dramatically following microbial infections, and this may be particularly helpful for the diagnosis and monitoring of bacterial septicemia in neonates and other immunocompromised patients at risk. In children, CRP is useful for differential diagnosis of bacterial and viral meningitis.

Because the biological half-life of this protein is only 24 hours, CRP accurately parallels the activity of the inflammation process and the CRP concentration decreases much faster than ESR or any other acute phase parameter, which is particularly useful in monitoring appropriate treatment of bacterial diseases with antibiotics.

C-Reactive Protein measurements during the early and late post transplant period of bone marrow and organ transplantations is particularly useful in the management of interfering infections in these immunosuppressed patients.

3. PRINCIPLE OF THE CRP ELISA

Microtiterstrips coated with anti-CRP antibody are incubated with diluted standard sera and patient samples. During this incubation step CRP is bound specifically to the wells. After removal of the unbound serum proteins by a washing procedure, the antigen-antibody complex in each well is detected with specific peroxidase-conjugated antibodies.

After removal of the unbound conjugate, the strips are incubated with a chromogen solution containing tetramethylbenzidin and hydrogen peroxide: a blue colour develops in proportion to the amount of immunocomplex bound to the wells of the strips. The enzymatic reaction is stopped by the addition of 0.5 M H₂SO₄ and the absorbance values at 450 nm are determined.

A standard curve is obtained by plotting the absorbance values versus the corresponding standard values. The concentration of CRP in patient samples is determined by interpolation from the standard curve.

4. REAGENTS

1. 1. Coated Microtiterstrips – MTP

12 x 8-well strips coated with monoclonal antibodies to human CRP.

2. Standard Serum – CAL N

5 vials, each containing 1/10 prediluted CRP standard solutions (0.2 ml), N having following values:

N=0: 0 µg/ml; N=5: 5 µg/ml; N=25: 25 µg/ml; N=50: 50 µg/ml;

N=100: 100 µg/ml. Contains 0.09 % NaN₃.

Calibrated against the NIBSC 1st International Standard, 85/506.

3. Conjugate – CONJ

1 vial, containing peroxidase conjugated monoclonal anti-human CRP antibodies (12 ml). Contains antimicrobial agents and an inert red dye.

4. Specimen Dilution Buffer – DIL 5x

1 vial, containing 40 ml dilution buffer 5x concentrated. Contains 0.09 % NaN₃ and an inert green dye.

5. Washing Solution – WASH 20x

1 vial containing 50 ml 20x concentrated phosphate buffered washing solution.

6. Chromogen Solution – CHROM

1 vial, containing 15 ml of a solution containing H₂O₂ and tetramethylbenzidin.

7. Stopping Solution – STOP

1 vial, containing 12 ml of 0.5 M H₂SO₄.

5. MATERIALS REQUIRED BUT NOT SUPPLIED


- Precision micropipettes and standard laboratory pipettes.
- Clean standard laboratory volumetric glassware.
- Clean glass or plastic tubes for the dilution of the samples.
- A microtiterplate reader capable of measuring absorbencies at 450 nm.

6. WARNINGS AND PRECAUTIONS FOR USERS

- For in vitro diagnostic use only.
- Human blood components used in the preparation of the standard sera have been tested and found to be nonreactive for hepatitis B surface antigen and HIV I. Since no known method can ever offer complete assurance that products derived from human blood will not transmit hepatitis or other viral infections, it is recommended to handle these standard sera in the same way as potentially infectious material. Dispose patient samples and all materials used to perform this test as if they contain infectious agents.
- Do not mix reagents or coated microtiterstrips from kits with different lot numbers.
- Some kit components contain sodium azide as a preservative. In order to prevent the formation of potentially explosive metal azides in laboratory plumbing, flush drains thoroughly after disposal of these solutions.

7. STORAGE CONDITIONS

- Store the microtiterstrips in their original package with the desiccant until all the strips have been used.
- Never use any kit components beyond the expiration date.

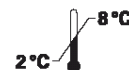

8. SPECIMEN COLLECTION AND PREPARATION

Human serum and plasma may be used in this assay. Remove serum from clot as soon as possible to avoid haemolysis. Lipemic and/or haemolysed samples can cause false results. Transfer the serum to a clean storage tube. Specimens may be stored at 2–8 °C for a few days, or they can be stored frozen for a longer period of time. Avoid repeated freezing and thawing.

9. ASSAY PROCEDURE

General Remarks:

- Use a separate disposable tip for each sample transfer to avoid cross-contamination.
- All reagents must be allowed to come to room temperature before use. All reagents must be mixed without foaming.
- Once the assay has been started, all steps should be completed without interruption.
- If an ELISA Washer is used, adaptation of the washing step might be necessary to obtain optimal results.

Reconstitution of the Reagents

Washing Solution:

Dilute 50 ml of concentrated Washing Solution (5) to 1000 ml with distilled water. Reconstituted solution can be stored at least 1 month, store at 2–8 °C.

At higher temperatures, the concentrated Washing Solution (5) may appear cloudy without affecting its performance. Upon dilution, the solution will be clear.

Specimen Dilution Buffer:

Dilute 40 ml of the concentrated Specimen Dilution Buffer (4) to 200 ml with distilled water. Reconstituted solution can be stored at least 3 months or as long as solution remains clear. Store at 2–8 °C.

Assay Procedure

1. The 10x prediluted standard sera (2) are diluted 1:100 as follows: pipette 10 µl of each calibrator into separate glass or plastic dilution tubes. Add 990 µl of diluted Specimen Dilution Buffer (4) and mix carefully.
2. The patient samples are diluted 1:1000 in two consecutive steps: pipette 10 µl of each patient sample into separate glass or plastic dilution tubes and add 990 µl of diluted Specimen Dilution Buffer (4). Mix thoroughly. Add 450 µl of diluted Specimen Dilution Buffer to 50 µl of these 100x prediluted samples. Mix thoroughly.

Warning: do not store the diluted samples for more than 8 hours.

3. Pipette 100 µl of the diluted calibrators and samples into each of a pair of adjacent wells (1).
4. Incubate the covered microtiterstrips for 30±2 min at room temperature.
5. Wash the microtiterstrips three times with Washing Solution. This can either be performed with a suitable microtiterplate washer or by briskly shaking out the contents of the strips and immersing them in washing solution. During the third step, the washing solution is left in the strips for 2–3 min. Change washing solution for each cycle. Finally empty the microtiterstrips and remove excess fluid by blotting the inverted strips on adsorbent paper.

6. Add 100 μ l of Conjugate Solution (3) and incubate the covered microtiterstrips for 30 ± 2 min at room temperature.
7. Repeat the washing procedure as described in 5.
8. Add 100 μ l of Chromogen Solution (6) to each well.
9. Incubate for 10 ± 2 min at room temperature. Avoid light exposure during this step.
10. Add 50 μ l of Stopping Solution (7) to each well.

Warning: determine the absorbance of each well at 450 nm within 30 min following the addition of acid.

10. CALCULATIONS

The average absorbance value of each calibrator is plotted against the corresponding CRP-value and the best calibration curve (e.g. log/linear) is constructed.

Use the average absorbance of each patient sample obtained in the CRP ELISA to determine the corresponding value by simple interpolation from the curve.

Depending on the experience and/or availability of computer capability, other methods of data reduction may be used.

The following specifications must be met for each run to be valid:

O.D. value for the zero calibrator: < 0.080

O.D. value for the highest value calibrator: > 1.000

If one of the specifications is not met, the test run should be repeated.

Example of typical O.D. values:

CALIBRATOR µg/ml	O.D. value
0	0.019
5	0.240
25	0.821
50	1.301
100	2.018

11. PERFORMANCE CHARACTERISTICS

- **Cross-reactivity**

The BioVendor CRP ELISA recognizes natural and recombinant human CRP. No cross-reactivity was observed with following factors, prepared at 1 µg/ml in sample diluent: human pentraxin 2; human pentraxin 3; human monomeric CRP; rat CRP.

- **Analytical Sensitivity**

The minimal detectable concentration is lower than 1 µg/ml.

- **Precision**

<i>Intra Assay (n=10)</i>	Level 1	Level 2		
Mean (µg/ml)	5.2	48.3		
SD (µg/ml)	0.27	3.3		
CV (%)	5.12	6.84		
<i>Inter Assay (n=7)</i>	Level 1	Level 2	Level 3	
Mean (µg/ml)	4.3	31.0	67.2	
SD (µg/ml)	0.6	3.6	8.5	
CV (%)	14.3	11.6	12.7	

12. CLINICAL DATA

Serum and plasma samples from 360 healthy blood donors were tested with the CRP ELISA, following distribution of CRP levels has been found:

Concentration Range CRP µg/ml	Number of donors
≤ 5	328
6–25	31
26–50	1
> 50	0

All individuals have small amounts of CRP in their blood. The upper limit of the normal range is situated between 5 and 8 µg/ml (343 donors in the population of 360 individuals had CRP levels < 8 µg/ml).


13. TROUBLESHOOTING


In case of high background signal, the washing was insufficient. Repeat the test with more vigorous washing (increased number of cycles, soak time).

14. REFERENCES

1. POWELL LJ. C-Reactive Protein - a review, *Am J Med Technol*, 1979, 45:138-142.
2. GEWURZ H., MOLD C., SIEGEL J., FIEDEL B. C-Reactive Protein and the Acute Phase Response, *Adv Intern Med*, 1982, 27:345-372.
3. HELGESON NGP., ADAMSON DM., PIKE RB., JAMES DS., NICODEMUS DS., LEE BA., MILLER GW. C-Reactive Protein, *Laboratory Medicine*, 1973, Vol. 2 (Race GJ., Ed.), Harper & Row, Hagerstown, chapter 29.
4. JOHNSON HL., CHIOU CC., CHO CT. Applications of acute phase reactants in infectious diseases, *J Microbiol Immunol Infect*, 1999, 32(2):73-82.

NOTES

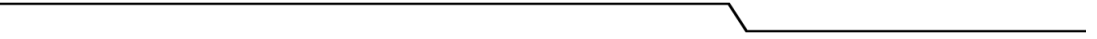

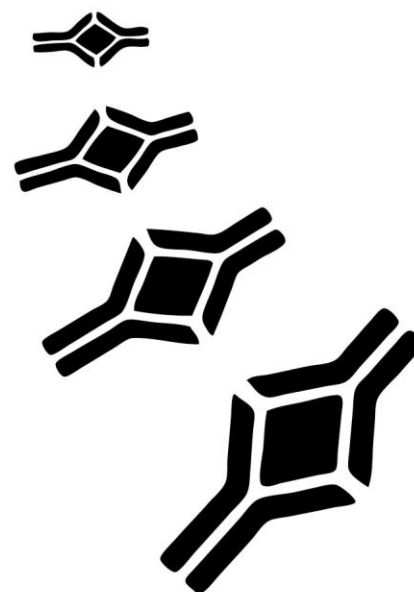

BioVendor – Laboratorní medicína, a.s.

Karasek 1767/1, 621 00 Brno, Czech Republic

Phone: +420-549-124-185, Fax: +420-549-211-460

E-mail: info@biovendor.com, sales@biovendor.com

Web: www.biovendor.com

There are BioVendor branches and distributors near you.

To find the office closest to you, visit www.biovendor.com/contact



Date of last revision: 5.12.2016