




Instruction for use


ADMA - ELISA

Enzyme Immunoassay
for the Quantitative Determination of
Endogenous Asymmetric Dimethylarginine (ADMA)
in Serum or Plasma



REF EA201/96

 12 x 8

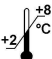 2 – 8 °C

DLD Gesellschaft für Diagnostika und medizinische Geräte mbH
Adlerhorst 15 • 22459 Hamburg • Germany
Tel +49-40-555 87 10 • Fax +49-40-555 87 111
Internet: <http://www.dld-diagnostika.de> • E-Mail: contact@dld-diagnostika.de

Contents

1. Introduction and Principle of the Test	Page	3
2. Precautions	Page	3
3. Storage and Stability	Page	4
4. Contents of the Kit	Page	4
5. Sample Collection	Page	5
6. Preparation of Reagents and Samples	Page	6
7. Test Procedure (ELISA)	Page	8
8. Calculation of the Results	Page	9
9. Assay Characteristics	Page	10
10. Literature	Page	13
Pipetting Scheme Sample Preparation	Page	15
Pipetting Scheme ELISA	Page	16

1. Introduction and Principle of the Test

The vascular endothelium plays a central role in the regulation of vascular structure and function, mainly due to the formation of endothelium-derived nitric oxide (NO). NO has been named an “endogenous anti-atherogenic molecule” due to its diverse regulatory functions in vascular homeostasis.

NO is formed by the enzyme NO synthetase (NOS) from the amino acid precursor L-arginine. NOS activity can be downregulated by asymmetric dimethylarginine (ADMA), an endogenous inhibitor of NOS.

The effects of ADMA on NO synthesis and NO-mediated pathophysiological processes have been described in numerous experimental studies. Moreover, elevated ADMA levels in plasma have been found in clinical studies including patients with hypercholesterolemia, hypertension, chronic heart failure, chronic renal failure and other internal disorders.

Recent prospective and cross-sectional studies indicated that elevated ADMA levels are a risk factor for future cardiovascular events and total mortality. ADMA may have diagnostic relevance as a novel cardiovascular risk marker.

The competitive ADMA-ELISA uses the microtiter plate format. ADMA is bound to the solid phase of the microtiter plate. ADMA in the samples is acylated and competes with solid phase bound ADMA for a fixed number of rabbit anti-ADMA antiserum binding sites. When the system is in equilibrium, free antigen and free antigen-antiserum complexes are removed by washing. The antibody bound to the solid phase ADMA is detected by anti-rabbit/oxidase. The substrate TMB / oxidase reaction is monitored at 450 nm. The amount of antibody bound to the solid phase ADMA is inversely proportional to the ADMA concentration of the sample.

2. Precautions

- For in vitro use only.
- Disposable gloves should be used.
- Material of animal origin used in the preparation of the kit has been obtained from animals certified as healthy but these materials should be handled as potentially infectious.

3. Storage and Stability

On arrival, store the kit at 2-8 °C. Once opened the kit is stable until its expiry date. For stability of prepared reagents refer to Preparation of Reagents.

Do not use components beyond the expiration date shown on the labels. Do not mix various lots of any kit component within an individual assay.

4. Contents of the Kit

4.1 **MT-Strips** **STRIPS** 12 strips
8 wells each, break apart
precoated with ADMA

4.2 **Standards 1 - 6** **CAL 1 - 6** 6 vials
Each 4 ml, ready for use
Concentrations:

Standard	1	2	3	4	5	6
µmol/l	0	0.1	0.3	0.6	1.0	5.0

4.3 **Control 1 & 2** **CON 1 & 2** 2 vials
Each 4 ml, ready for use
Range: see q.c. certificate

4.4 **Acylation Buffer** **ACYL-BUFF** 1 bottle
3.5 ml, ready for use

4.5 **Acylation Reagent** **ACYL-REAG** 2 vials
lyophilised, dissolve content
in 2.8 ml Solvent before use; if required
combine the contents of both vials (see also 6.4.)

4.6 **Antiserum** **AS** 1 vials
5.5 ml, ready for use
Rabbit-anti-N-acyl-ADMA

4.7 **Enzyme Conjugate** **CONJ** 1 vial
12 ml, ready for use
goat anti-rabbit-IgG-peroxidase

4.8 **Wash Buffer** **WASH** 1 bottle
20 ml, concentrated
Dilute content with dist. water to 500 ml total volume.

4.9	Substrate 12 ml TMB solution, ready for use	SUB	1 vial
4.10	Stop Solution 12 ml, ready for use Contains 0.3 M sulphuric acid, not corrosive	STOP	1 vial
4.11	Reaction Plate for acylation	ACYL-PLATE	1 piece
4.12	Equalizing Reagent lyophilized, dissolve content with 20.5 ml dist. water, dissolve carefully to minimize foam formation	EQUA-REAG	1 vial
4.13	Solvent 6 ml, contains acetone/ DMSO (please note that Solvent reacts with many plastic materials including plastic trays; Solvent does not react with normal pipette tips and with glass devices)	SOLVENT	1 vial

Additional materials and equipment required but not provided:

- Pipettes (20, 25, 50, 100 and 200 µl)
- Orbital shaker
- Microplate washing device
- Microplate photometer (450 nm)
- Vortex mixer
- Roll mixer

5. Sample Collection

5.1. Serum and Plasma

The test can be performed with serum as well as with EDTA plasma.

Hemolytic and lipemic samples should not be used.

The samples can be stored up to 24 hours at 2 - 8 °C. For a longer storage (up to 24 months) the samples must be frozen at -20 °C

Repeated freezing and thawing should be avoided.

6. Preparation of Reagents and Samples

6.1. Microtiter strips **STRIPS**

Before opening the packet of strip wells, allow it to stand at room temperature for at least 10 minutes. After opening, keep any unused wells in the original foil packet with the desiccant provided. Reseal carefully and store at 2-8 °C.

6.2 Wash Buffer **WASH**

Dilute the content with dist. water to a total volume of 500 ml. The diluted wash buffer has to be stored at 2 - 8 °C.

6.3. Equalizing Reagent **EQUA-REAG**

Dissolve the content with 20.5 ml dist. water, mix shortly and leave on a roll mixer or orbital shaker for 30 minutes. Handle carefully in order to minimize foam formation. The reconstituted Equalizing Reagent should be stored frozen at -20 °C and is stable for a minimum of 1 year.

6.4. Acylation Reagent **ACYL-REAG**

Dissolve the content of one bottle in 2.8 ml Solvent and shake for 5 minutes on a orbital shaker. After use the reagent has to be discarded. The Acylation Reagent has always to be prepared immediately before use. The second bottle allows a second run of the test. If the whole kit is to be used in one run it is recommended to pool the dissolved contents of the two vials of Acylation Reagent.

Please note that Solvent reacts with many plastic materials including plastic trays. Solvent does not react with normal pipette tips and with glass devices.

Attention

Solvent is volatile and the dissolved Acylation Reagent evaporates quickly. Therefore, please do not use a tray with big surface together with a multichannel pipette for pipetting Acylation Reagent. Rather, use an Eppendorf multipette, or similar device, fill the syringe directly from the vial (using a yellow tip) with dissolved Acylation Reagent and add well by well.

All other reagents are ready for use.

6.5. Preparation of Samples (Acylation)

The wells of the reaction plate for the acylation can be used only once. Please mark the respective wells before use to avoid repeated use.

1. Pipette each 20 μ l standard 1 - 6, each 20 μ l control 1 & 2 and each 20 μ l patient sample into the respective wells of the Reaction Plate.
2. Pipette 25 μ l Acylation Buffer into all wells.
3. Pipette 200 μ l Equalizing Reagent into all wells.
4. Mix the reaction plate for 10 seconds.
5. Prepare Acylation Reagent just before use and pipette 50 μ l prepared Acylation Reagent each into all wells, mix **immediately**.

Attention

Solvent is volatile and the dissolved Acylation Reagent evaporates quickly. Therefore, please do not use a tray with big surface together with a multichannel pipette for pipetting Acylation Reagent. Rather, use an Eppendorf multipipette with a yellow tip (or similar device), fill the syringe directly from the vial with dissolved Acylation Reagent and add well by well.

6. Incubate for 90 minutes at room temperature (approx. 20 $^{\circ}$ C) on an orbital shaker. Do not cover the wells or the plate; leave the plate open on the shaker.

Take each 50 μ l for the ADMA-ELISA.

7. Test Procedure ELISA

Bring all reagents to room temperature and mix them carefully, avoid development of foam.

7.1 Sample Incubation

Pipette each 50 μ l prepared Standards 1 to 6, 50 μ l prepared controls and 50 μ l prepared samples into the respective wells of the coated microtiter strips (duplicates are recommended).

Pipette each 50 μ l Antiserum into all wells and shake shortly on an orbital shaker.

Cover the plate with adhesive foil and incubate Microtiter Strips for 15 –20 hours (overnight) at 2 – 8 $^{\circ}$ C.

7.2 Washing

Discard or aspirate the contents of the wells and wash thoroughly with each 250 μ l Wash Buffer (Shake shortly on an orbital shaker). Repeat the washing procedure 4 times. Remove residual liquid by tapping the inverted plate on clean absorbent paper.

7.3 Conjugate Incubation

Pipette each 100 μ l enzyme conjugate into all wells.

Incubate for 60 minutes at room temperature on an orbital shaker.

7.4 Washing

Repeat step 7.2.

7.5 Substrate Incubation

Pipette each 100 μ l Substrate into all wells and incubate for 20 to 30 minutes at room temperature on an orbital shaker.

7.6 Stopping

Pipette each 100 μ l Stop Solution into all wells.


7.7 Reading

Read the optical density at 450 nm (reference wavelength between 570 and 650 nm) in a microplate photometer.

8. Calculation of the Results

On a semilogarithmic graph paper the concentration of the standards (x-axis, logarithmic) are plotted against their corresponding optical density (y-axis, linear). Cubic spline, 4 parameter or similar iteration procedures are recommended for evaluation of the standard curve. The concentration of the controls and samples can be read directly from this standard curve by using their average optical density.

Typical standard curve:

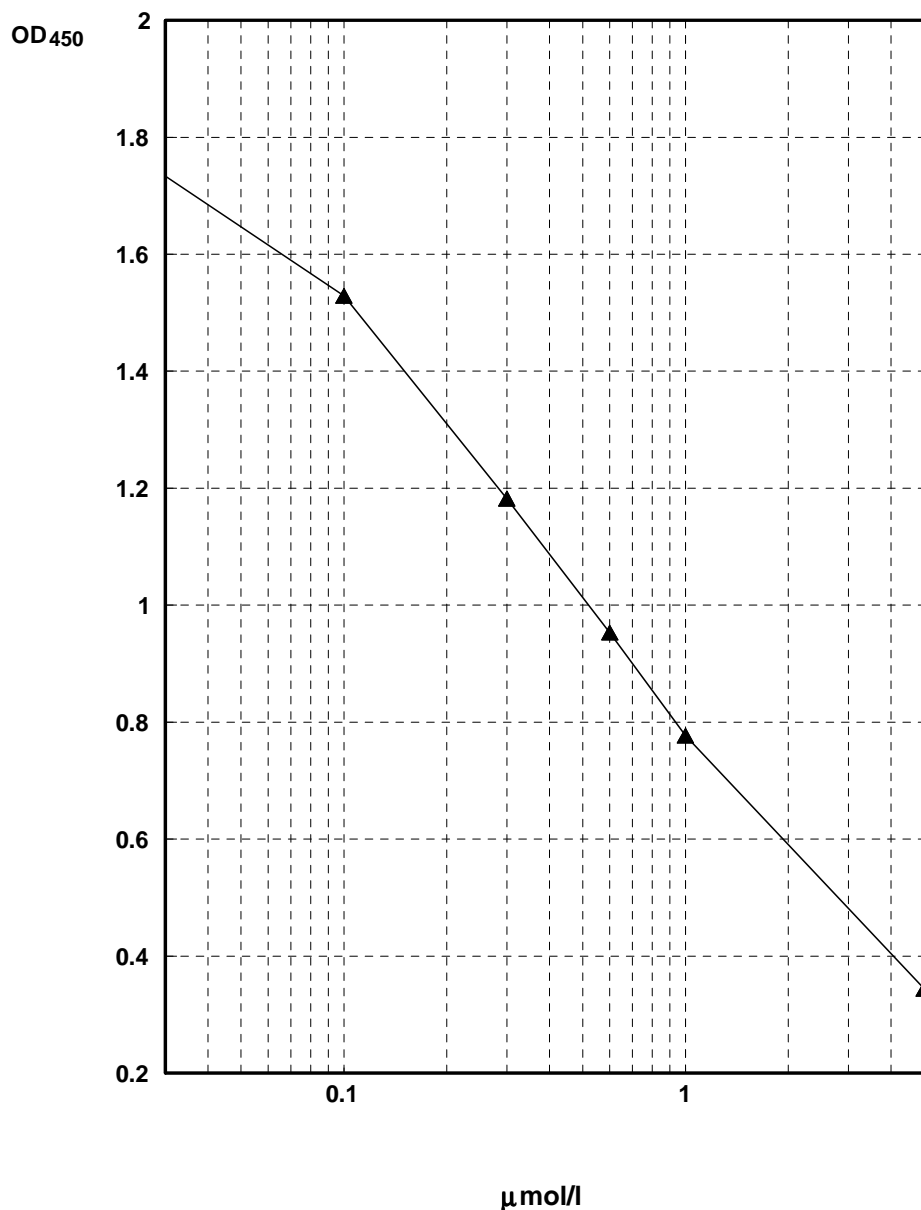

9. Assay Characteristics

Expected Values


0.4 – 0.75 µmol/l (80 – 150 ng/ml)

The reference ranges given above should only be taken as a guideline. It is recommended that each laboratory should establish its own reference values.

Sensitivity

0.05 µmol/l

Correlation to LC-MS-MS

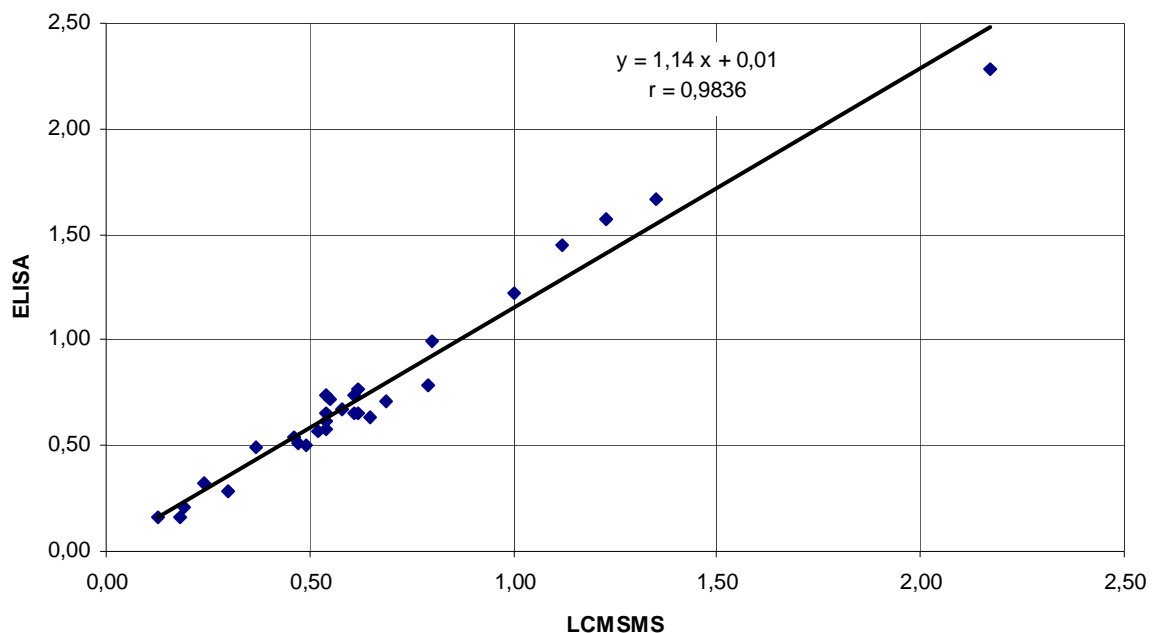

The figure shows the correlation to the LC-MS-MS method. More data with higher number of patient samples are evaluated at the moment. Comparative investigations to HPLC without subsequent MS lead to inconsistent results.

Specificity (Cross Reactivity)

Structural related components were tested for possible interference with the antisera against ADMA used in the ELISA method. The tested compounds were Arginine, Monomethylarginine (NMMA) und SDMA.

Substance	ED-50-Value (ng/ml)	Cross Reactivity (%)
ADMA	126	100
Arginin	660,000	< 0.02
NMMA	12,200	1.0
SDMA	10,500	1.2

Recovery

Increasing amounts of ADMA were added to a serum sample. Each spiked sample was assayed. The analytical recovery of ADMA was estimated at twelve different concentrations by using the theoretically expected and the actually measured values. The mean recovery from all concentrations was 97% (90 - 104%).

added	measured	expected	% recovery
0	0.43		
0.1	0.55	0.53	104
0.19	0.61	0.62	98
0.28	0.67	0.71	94
0.37	0.79	0.80	99
0.45	0.83	0.88	97
0.73	1.08	1.16	93
0.97	1.33	1.40	95
1.20	1.60	1.63	98
1.42	1.67	1.85	90
1.92	2.24	2.35	95
2.33	2.76	2.76	100
2.80	3.32	3.23	103

mean value 97

concentrations in $\mu\text{mol/l}$

Linearity

The linearity of the ELISA method was investigated using nine different dilutions of a serum sample. The mean linearity from all dilutions was 97% (88 - 112%)

dilution	measured	recalculated value	recovery %
orig.	2.819		
3 + 1	2.090	2.787	99
2 + 1	1.745	2.618	93
1 + 1	1.236	2.472	88
1 + 2	0.860	2.580	91
1 + 4	0.535	2.675	95
1 + 7	0.361	2.888	102
1 + 10	0.286	3.146	112
1 + 15	0.178	2.848	101
1 + 25	0.104	2.704	96

mean recovery

97

concentrations in $\mu\text{mol/l}$

Reproducibility

The reproducibility of the ELISA method was investigated by determining the intra- und inter-assay-coefficients of variation (cv) by repeated measurements of different serum samples with different ADMA concentrations.

Intra-Assay Variation

sample	n =	mean value	sd	cv (%)
1	42	0.66	0.037	5.7
2	42	1.01	0.066	6.4

Inter-Assay Variation

sample	n =	mean value	sd	cv (%)
1	28	0.63	0.07	10.3
2	28	1.01	0.10	9.8
3	28	1.38	0.13	9.4
4	28	2.26	0.19	8.3

9. Literature

Literature using the ADMA-ELISA from DLD Diagnostika

Schulze F, Wesemann R, Schwedhelm E, Sydow K, Albsmeier J, Cooke JP, Böger RH.

Determination of ADMA using a novel ELISA assay.

Clin. Chem. Lab. Med. 2004; 42: 1377-1383

Krempf TK, Kähler J, Maas R, Silberhorn L, Meinertz T, Böger RH.

Elevation of asymmetric dimethylarginine (ADMA) in patients with unstable angina and recurrent cardiovascular events.

Eur. Heart J. 2005; 26: 1846-1851

Schulze F, Maas R, Freese R, Schwedhelm E, Silberhorn L, Böger RH.

Determination of a reference value for N,N-dimethyl-L-arginine in 500 subjects.

Eur. J. Clin. Invest. 2005; 35 : 622-626

Schnabel R, Blankenberg S, Lubos E, Lackner KJ, Rupprecht HJ, Espinola-Klein C, Jachmann N, Post F, Peetz D, Bickel C, Cambien F, Tiret L, Münzel T.

Asymmetric dimethylarginine and the risk of cardiovascular events and death in patients with coronary artery disease: results from the AtheroGene Study.

Circ. Res. 2005; 97: e53-59

O'Dwyer MJ, Dempsey F, Crowley V, Kelleher D, McManus R, Ryan T.

Septic shock correlates with ADMA levels which may be influenced by a polymorphism in DDAH II: a prospective observational study.

Crit. Care 2006; 10: (5): R139

Antoniades C, Tousoulis D, Marinou K, Vasiliadou C, Tentolouris C, Bouras G, Pitsavos C, Stefanidis C.

Asymmetrical dimethylarginine regulates endothelial function in methionine-induced but not in chronic homocystinemia in humans: effect of oxidative stress and proinflammatory cytokines

Am. J. Clin. Nutr. 2006; 84: 781-788

Wang TZ., Chen WJ., Cheng WC., Lin JW., Chen MF., Lee YT.

Relation of improvement in endothelium-dependent flowmediated vasodilation after Rosiglitazone to changes in asymmetric dimethylarginine, endothelin-1, and C-reactive protein in nondiabetic patients with the metabolic syndrome

Am. J. Cardiol. 2006; 9: 1057-1062

Wanby P., Nilsson I., Brudin L., Nyhammar I., Gustafsson I., Carlsson M.

Increased plasma levels of asymmetric dimethylarginine in patients with carotid stenosis: no evidence for the role of the common FABBP2 A54T gene polymorphism

Acta Neurol. Scand. 2007; 115: 90-96

Konishi H, Sydow K, Cooke JP.

Dimethylarginine dimethylaminohydrolase promotes endothelial repair after vascular injury

J. Am. Coll. Cardiol. 2007; 49: 1099-1105

Iribarren C, Husson G, Sydow K, Wang BY, Sidney S, Cooke JP.

Asymmetric dimethyl-arginine and coronary artery calcification in young adults entering middle age: the CARDIA Study

Eur. J. Cardiovasc. Prev. Rehabil. 2007; 14:222-229

Melikian N, Wheatcroft SB, Ogah OS, Murphy C, Chowienczyk PJ, Wierzbicki AS, Sanders TA, Jiang B, Duncan ER, Shah AM, Kearney MT.

Asymmetric dimethylarginine and reduced nitric oxide bioavailability in young Black African men

Hypertension 2007; 49: 873-877

Horowitz JD, Heresztyn T.

An overview of plasma concentrations of asymmetric dimethylarginine (ADMA) in health and disease and in clinical studies: Methodological considerations.

J. Chromatogr. B Analyt. Technol. Biomed. Life Sci. 2007; epub ahead of print

Korish AA, Arafah MM.

Catechin combined with vitamins C and E ameliorates insulin resistance (IR) and atherosclerotic changes in aged rats with chronic renal failure (CRF)

Arch. Gerontol. Geriatr. 2007; in press

Charitidou C, Farmakiotis D, Zournatzi V, Pidonia I, Pegiou T, Karamanis N, Hatzistilianou M, Katsikis I, Panidis D.

The administration of estrogens, combined with anti-androgens, has beneficial effects on the hormonal features and asymmetric dimethyl-arginine levels, in women with the polycystic ovary syndrome

Atherosclerosis 2007; in press

General Literature

Vallance P, Leone A, Calver A, Collier J, Moncada S.

Accumulation of an endogenous inhibitor of NO synthesis in chronic renal failure

Lancet 1992; 339: 572 - 575

Stühlinger M, Abbasi F, Chu JW, Lamendola C, McLaughlin TL, Cooke JP, Reaven GM, Tsao PS.

Relationship between insulin resistance and an endogenous nitric oxide synthase inhibitor

J. Am. Med. Assoc. 2002; 287: 1420-1426

Zoccali C, Bode-Böger SM, Mallamaci F, Benedetto FA, Tripepi G, Malatino L, Cataliotti A, Bellanuova I, Fermo I, Frölich JC, Böger RH.

Asymmetric dimethylarginine (ADMA): An endogenous inhibitor of nitric oxide synthase predicts mortality in end-stage renal disease (ESRD)

Lancet 2001; 358: 2113-2117

Nijveldt RJ, Teerlink T, Van der Hoven B, Siroen MP, Kuik DJ, Rauwerda JA, van Leeuwen PA.

Asymmetrical dimethylarginine (ADMA) in critically ill patients: high plasma ADMA concentration is an independent risk factor of ICU mortality

Clin. Nutr. 2003; 22: 23-30

Savvidou MD, Hingorani AD, Tsikas D, Frolich JC, Vallance P, Nicolaidis KH.

Endothelial dysfunction and raised plasma concentrations of asymmetric dimethylarginine in pregnant women who subsequently develop pre-eclampsia

Lancet 2003; 361: 1511-1517

Böger RH.

The emerging role of asymmetric dimethylarginine as a novel cardiovascular risk factor

Cardiovasc. Res. 2003; 59: 824-833

Lu TM, Ding YA, Lin SJ, Lee WS, Tai HC.

Plasma levels of asymmetrical dimethylarginine and adverse cardiovascular events after percutaneous coronary intervention.

Eur Heart J. 2003; 24: 1912-1919

Pipetting Scheme Sample Preparation

		Standards	Control	Sample
Standard 1 - 6	μl	20		
Control 1 & 2	μl		20	
Patient Sample	μl			20
Acylation Buffer	μl	25	25	25
Equalizing Reagent	μl	200	200	200

shake for 10 seconds

freshly prepared Acylation Reagent	μl	50	50	50
---------------------------------------	----	----	----	----

incubate for 90 minutes at room temperature on an orbital shaker
do not cover wells or plate, leave the plate open on the shaker

Pipetting Scheme ELISA

		Standard	Control	Patient-Sample
Standard 1 - 6	µl	50		
Control 1 & 2	µl		50	
Patient Sample	µl			50
Antiserum	µl	50	50	50

shake shortly on an orbital shaker

incubate 15 – 20 hours (overnight) at 2 - 8 °C covered with foil

wash 4 x with each 250 µl Wash Buffer

Enzyme Conjugat	µl	100	100	100
-----------------	----	-----	-----	-----

shake for 60 minutes at room temperature

wash 4 x with each 250 µl Wash Buffer

Substrate	µl	100	100	100
-----------	----	-----	-----	-----

shake for 20 - 30 minutes at room temperature

Stop Solution	µl	100	100	100
---------------	----	-----	-----	-----

read absorbance at 450 nm