

DRG[®] Homocysteine (EIA-2925)



Revised 27 Feb. 2014 rm (Vers. 5.1)

USA: RUO

This kit is intended for Research Use Only.

Not intended for use in diagnostic procedures.

Intended Use

The Homocysteine ELISA is intended for measurement of total L-homocysteine in human serum or plasma.

Assay Principle

The Homocysteine ELISA is an enzyme immunoassay for the determination of Hcy in blood.[16] Protein-bound Hcy is reduced to free Hcy and enzymatically converted to S-adenosyl-L-homocysteine (SAH) in a separate procedure prior to the immunoassay.[17] The enzyme is specific for the L-form of homocysteine, which is the only form present in the blood.


Reduction

Hcy, mixed disulfide and protein-bound forms of Hcy in the sample are reduced to free Hcy by use of dithiothreitol (DTT).



Enzymatic conversion

Hcy in the test sample is converted to S-adenosyl-L-homocysteine by the use of SAH hydrolase and excess adenosine (Ad).

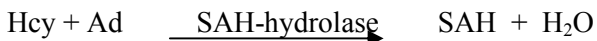

The following solid-phase enzyme immunoassay is based on competition between SAH in the sample and immobilised SAH bound to the walls of the microtitre plate for binding sites on a monoclonal anti-SAH antibody. After removal of unbound anti-SAH antibody, a secondary rabbit anti-mouse antibody labelled with the enzyme horse radish peroxidase (HRP) is added. The peroxidase activity is measured spectrophotometrically after addition of substrate, and the absorbance is inversely related to the concentration of Hcy in the sample.

Warnings and Precautions

1. Reagent D contains 0.15% merthiolate ($\leq 0.074\%$ mercury), and is classified as “Harmful”. Please handle and dispose of properly (See section “Product Safety Information”).
2. 0.01% merthiolate is used as preservative in some reagents. Each kit contains less than 0.028% mercury. Please handle and dispose of appropriately.
3. Reagent F contains mouse antibodies and Reagent G contains rabbit antibodies.

DRG® Homocysteine (EIA-2925)



Revised 27 Feb. 2014 rm (Vers. 5.1)

USA: **RUO**

4. Reagent S contains 0.8M sulphuric acid, and is classified as “Irritant”. Please handle and dispose of properly (See section “Product Safety Information”).
5. Calibrators, Controls, Reagent A and Reagent E contain less than 0.10% sodium azide as preservative. Sodium azide may react with lead and copper plumbing to form highly explosive metal azides. On disposal, flush with a large volume of water to prevent azide build-up.
6. Controls contain sera originating from human blood samples. The source materials have been tested and found to be negative for Hepatitis B Surface Antigen (HBsAG), HIV-1 Antigen (HIVAg), HCV antibody, HIV-1/2 antibody, HTLV-1/2 antibody and Hepatitis B core Antibody (HBc). However, blood derivatives should be handled according to recommended procedures for handling infectious material. HHS publication no. (CDC) 93-8395 [18] or local/national guidelines on laboratory safety procedures should be consulted.
7. Reagents with different lot numbers must not be interchanged.
8. Do not use the kit after the expiration date on the outer box.

Kit Components**Homocysteine ELISA Kit, 96 wells**

Symbol	Kit Components	Solution	Component Description	Volume
BUF	[REAG A]	Assay Buffer	Phosphate buffer, sodium azide	54 mL
ADENO DTT	[REAG B]	Adenosine/DTT	Adenosine / dithiothreitol, citric acid	3.5 mL
SAH HYDROL	[REAG C]	SAH-hydrolase	Recombinant S-adenosyl-L-homocysteine hydrolase, trisbuffer, glycerol, methylparaben	3.5 mL
ENZ INH	[REAG D]	Enzyme Inhibitor	Merthiolate, phosphate buffer	55 mL
ADENO DEAM	[REAG E]	Adenosine deaminase	Adenosine deaminase, phosphate buffer, sodium azide, BSA, phenol-red dye	55 mL
SAH Ab	[REAG F]	a-SAH Antibody	Monoclonal mouse-anti-S-adenosyl-L-homocysteine antibody, BSA, merthiolate	25 mL
ENZ CONJ	[REAG G]	Enzyme Conjugate	Rabbit anti-mouse-antibody enzyme conjugate, BSA, horse radish peroxidase, blue dye	15 mL
SUB TMB	[REAG H]	Substrate Solution	N-methyl-2-pyrrolidon, propyleneglycol	15 mL
STOP SOLN	[REAG S]	Stop Solution	0.8 M Sulphuric acid	20 mL
BUF WASH 10x	[BUF WASH]	Wash Buffer	Phosphate buffer, merthiolate, Tween 20, BSA	60 mL
CAL 1 - CAL 6	[CAL 1] – [CAL 6]	Standards	S-adenosyl-L-homocysteine (2, 4, 8, 15, 30, 50 µmol/L) in buffer with preservative	6 x 1.5 mL
SORB MT	[MICROTITER STRIPS]	Microtiter wells	Coated with S-adenosyl-L-homocysteine	12 x 8 wells

DRG[®] Homocysteine (EIA-2925)

Revised 27 Feb. 2014 rm (Vers. 5.1)

USA: **RUO**

Wash Buffer is concentrated and must be diluted (1+9) with purified water before use.
All other components are ready to use.

Homocysteine ELISA CONTROL kit (EIA-3329)

Kit Components	Component Description	Volume
[CONTROL L]	7.0 µmol/L homocysteine in diluted serum samples of human origin, phosphate buffer and preservative	1.5 mL
[CONTROL M]	12.5 µmol/L homocysteine in diluted serum samples of human origin, phosphate buffer and preservative	1.5 mL
[CONTROL H]	25.0 µmol/L homocysteine in diluted serum samples of human origin, phosphate buffer and preservative	1.5 mL

All controls are ready to use.

Homocysteine ELISA Wash Buffer

Kit Components	Component Description	Volume
[BUF WASH]	Phosphate buffer, merthiolate, Tween 20, BSA	1000 mL

[BUF WASH] is concentrated and must be diluted (1+9) with purified water before use.

Materials required but not provided in the kit:

- Homocysteine controls (see section “Quality Controls” for more information)
- Plastic or glass tubes for pre-treatment of samples
- Pipettes / multipipettes 25 µL, 100 µL, 200 µL and 500 µL or 8 channel multipipette for 100 µL and 200 µL
- Volumetric flask 50 mL and 600 mL
- Incubator, 37 °C
- Washer and reader (450 nm) for microtiter plates

Preparation and Storage of Kit Components

1. Components should be stored refrigerated (2 °C - 8 °C). Store all bottles upright and tightly capped. The components are stable until the stated expiration date when stored and handled as directed. Once the components in the Homocysteine ELISA Kit are opened they are stable for 12 weeks when stored at 2 °C - 8 °C.

DRG[®] Homocysteine (EIA-2925)**Revised 27 Feb. 2014 rm (Vers. 5.1)****USA: RUO**

2. The sample pre-treatment solution has to be made by mixing Reagent A, B and C (see Section "Procedure"). The solution is stable for one hour and has to be freshly made for each run.
3. The Wash Buffer must be diluted (1+9) with distilled water before use. The prepared Wash buffer is stable for 4 weeks when stored at room temperature (18 °C - 25 °C).
4. Reagent D and H are stored in dark bottles to avoid exposure to light.
5. It is important that the microtiter strips are kept dry, i.e. in the sealed bag with drying capsules, and stored refrigerated. Equilibration for a minimum of two hours is required to reach room temperature (18 °C - 25 °C). Leave the strips in the bag during equilibration.
6. Only the necessary number of microtiter strips should be kept in the frame during the run. Unused strips should be kept in the sealed bag with drying capsules.
7. Avoid exposure of the kit to temperatures exceeding 37 °C as this may denature the enzymes.

Specimen Collection and Preparation

EDTA-plasma or serum may be used with the Homocysteine ELISA assay.

As synthesis of Hcy continues in red blood cells after drawing, it is very important to prepare specimens as follows:

- Serum samples should be allowed to clot for no more than 30 minutes before centrifugation and separation of serum. Serum samples should be kept on ice prior to separation.
- EDTA-plasma samples must be centrifuged or put on ice immediately after drawing. EDTA-plasma samples may be kept on ice for up to 6 hours prior to separation by centrifugation.

Food consumption can affect circulating homocysteine levels. Protein rich meals give higher homocysteine values and should be avoided late in the day before sampling.[19, 20]

Standardised sampling procedures are crucial due to the above mentioned influencing factors. Complete mixing of thawed samples is required before use.

Plasma or serum samples may be stored for 12 weeks at 2 °C - 8 °C, for up to 3 weeks at room temperature (18 - 25 °C) and have been shown to be stable for at least 8 months if frozen at minus 20 °C.

Limitations

- If an automatic pipetting station is used, thorough washing of the tubing after addition of the blue coloured Reagent G may be needed - preferably with diluted acid followed by water. Any remaining solution in the tubing will interfere with the next assay step; i.e. addition of Reagent H.
- The washing procedure is critical for obtaining good precision. If manual washing is required, use 4 times 350 µL instead of 3 times 400 µL. After washing, empty the wells on paper towels.
- Avoid exposure of the kit to temperatures exceeding 37 °C as this may denature the enzymes.
- Specimens from individuals who are on drug therapy involving S-adenosyl-methionine may show falsely elevated levels of homocysteine.
- Specimens from individuals who have received preparations of mouse monoclonal antibodies may contain human anti-mouse antibody (HAMA). HAMA, present in serum or plasma specimens, may interfere in immunoassays which utilise mouse monoclonal antibodies. These specimens should not be assayed with the Homocysteine ELISA assay.

DRG[®] Homocysteine (EIA-2925)

Revised 27 Feb. 2014 rm (Vers. 5.1)

USA: **RUO**

- Specimens from individuals taking methotrexate, carbamazepine, phenytoin, nitrous oxide, anti-convulsants or 6-azauridine triacetate, may have elevated levels of homocysteine due to metabolic interference with the homocysteine metabolism.

Procedure

Make sure all solutions and microtiter strips are equilibrated to room temperature before use. Leaving the kit at room temperature over night is recommended. We recommend running the standards in duplicate and to performing a new calibration curve for each run to avoid run-to-run variations using coated microtiter plates.

Sample pre-treatment procedure

1. Sample pre-treatment solution must be made up no more than 1 hour prior to the start of the assay. Volume needed per 10 samples (no dead volume calculated):
4.5 mL REAG A
0.25 mL REAG B
0.25 mL REAG C Mix.
2. Dilute calibrators and samples/calibrators in plastic or glass tubes as follows:
25 µL standard/sample/calibrator
+ 500 µL sample pre-treatment solution Mix well.
Incubate for 30 minutes at 37 °C (Cap the tubes or cover with parafilm during incubation).
Note: Proceed with step 3 before the samples have cooled.
3. Add 500 µL REAG D
Mix well. Incubate for 15 minutes at 18 °C - 25 °C.
4. Add 500 µL REAG E
Mix well. Incubate for 5 minutes at 18 °C - 25 °C.

DRG[®] Homocysteine (EIA-2925)

Revised 27 Feb. 2014 rm (Vers. 5.1)

USA: RUO***Microtitre plate procedure***

5. Pipette 25 µL diluted calibrator / sample / calibrator from step 4 into the wells of the SAH-coated microtitre strips.
6. Add 200 µL REAG F to each well. Incubate for 30 min at 18 °C - 25 °C.
Use the enclosed lid during all incubations.
7. Wash with diluted Wash buffer.
Use 3 x 400 µL. If manual washing is required, use 4 times 350 µL instead of 3 times 400 µL. After washing, empty the wells on paper towels.
8. Add 100 µL REAG G to each well. Incubate for 20 min at 18 °C - 25 °C.
9. Wash with diluted Wash buffer.
Use 3 x 400 µL. If manual washing is required, use 4 times 350 µL instead of 3 times 400 µL. After washing, empty the wells on paper towels.
10. Add 100 µL REAG H to each well. Incubate for 10 min at 18 °C - 25 °C.
11. Add 100 µL REAG S to each well.
12. Shake and read at 450 nm within 15 minutes (Automatic plate shaker is preferred to ensure proper mixing).

Interpretation of Results

We recommend that a four parameter logistic curve fit is used for preparing the calibration curve and calculation of unknown samples.

Quality Control

We recommend that each laboratory use a homocysteine calibrator with known value.

Reference Range

The reference range should be determined by each laboratory to confirm the characteristics of the population being tested. As a point of reference, the following data may be used until the laboratory has analysed a sufficient number of samples to determine its own reference range.

The Hcy concentration in plasma or serum of healthy individuals varies with age, gender, geographical area and genetic factors.

DRG[®] Homocysteine (EIA-2925)



Revised 27 Feb. 2014 rm (Vers. 5.1)

USA: **RUO**

Measuring Range

The range is from 2 to 50 µmol/L.

References

1. Malinow, MR. Plasma homocyst(e)ine and arterial occlusive diseases: A mini-review. *Clin Chem* 1995;41:173-176.
2. Ueland PM. Homocysteine species as components of plasma redox thiol status. *Clin Chem* 1995;41:340-342.
3. Perry IJ, Refsum H, Morris RW, Ebrahim SB, Ueland PM, Shaper AG. Prospective study of serum total homocysteine concentration and risk of stroke in middle-aged British men. *The Lancet* 1995;346:1395-1398.
4. Finkelstein JD. Methionine metabolism in mammals. *J Nutr Biochem* 1990;1:228-237.
5. Mudd SH, Levy HL, Skovby F. Disorder of transsulfuration. In: C.R. Scriver, A.L. Beaudet, W.S. Sly, et al. eds. *The metabolic basis of inherited disease* New York: McGraw-Hill; 1995:1297-1327.
6. Clarke R, Daly L, Robinson K, Naughten E, Cahalane S, Fowler B, Graham I. Hyperhomocysteinemia: An independent risk factor for vascular disease. *N Engl J Med* 1991;324:1149-1155.
7. Deloughery TG, Evans A, Sadeghi A, McWilliams J, Henner WD, Taylor LM Jr, Press RD. Common mutation in methylenetetrahydrofolate reductase. Correlation with homocysteine metabolism and late-onset vascular disease. *Circulation* 1996;94:3074-3078.
8. Schmitz C, Lindpaintner K, Verhoef P, Gaziano JM, Buring J. Genetic polymorphism of methylenetetrahydrofolate reductase and myocardial infarction. *Circulation* 1996;94:1812-1814.
9. Boushey CJ, Beresford SAA, Omenn GS, Motulsky AG. A Quantitative Assessment of Plasma Homocysteine as a Risk Factor for Vascular Disease, *JAMA* 1995;274:1049-1 057.
10. Allen, RH, Stabler SP, Savage DG, Lindenbaum J. Diagnosis of cobalamin deficiency I: Usefulness of serum methylmalonic acid and total homocysteine concentrations, *Am J Hematol* 1 990;34:90-98
11. Stabler, SP, Lindenbaum J, Allen RH. The Use of homocysteine and other metabolites in the specific diagnosis of vitamin B-12 deficiency, *J Nutr* 1996;126: 1266S-1 272S.
12. Guttormsen AB, Svarstad E, Ueland PM, Refsum H. Elimination of homocysteine from plasma in subjects with end-stage renal failure. *Irish J Med Sci* 1 995;1 64 (Suppl. 1 5):8-9
13. Bostom AG, Lathrop L. Hyperhomocysteinemia in end-stage renal disease: Prevalence, etiology, and potential relationship to arteriosclerotic outcomes. *Kidney Int* 1997;52:10-20.
14. Desouza C, Keebler M, McNamara DB, Fonseca V. Drugs affecting homocysteine metabolism. Review. *Drugs* 2002;62 (4):605-616.
15. Refsum H, Ueland PM. Clinical significance of pharmacological modulation of homocysteine metabolism. Review. *Trends Pharmacol Sci* 1990; 11:411-416
16. Frantzen F, Faaren AL, Alheim I, Nordhei AK. An enzyme conversion immunoassay for determining total homocysteine in plasma or serum. *Clin Chem* 1998;44:31 1-316.
17. Sundrehagen E. Axis Biochemicals ASA. Enzymatic assay for homocysteine and a kit therefor. EP 623174/US5631 127.

DRG[®] Homocysteine (EIA-2925)

Revised 27 Feb. 2014 rm (Vers. 5.1)

USA: **RUO**

18. US Department of Health and Human Services. Biosafety in Microbiological and Biomedical Laboratories. HHS Publication No. (CDC) 93-8395. Washington, DC: US Government Printing Office, May 1999 (4th edition).
19. Ubbink JB, Vermaak Hayward WJ, van der Merwe A, Becker PJ. The effect of blood sample aging and food consumption on plasma total homocysteine levels. *Clinica Chimica Acta* 1992;207:1 19-1 28.
20. Guttormsen AB, Schneede J, Fiskerstrand T, Ueland PM, Refsum H. Plasma concentrations of homocysteine and other aminothiols compounds are related to food intake in healthy human subjects. *J Nutr* 1994;124:1934-1941.
21. Ueland PM, Refsum H, Stabler SP, Malinow R, Andersson A, Allen RH. Total homocysteine in plasma or serum: Methods and clinical applications. *Clin Chem* 1993;39:1764-1779.
22. Nehler MR, Taylor LM Jr, Porter JM. Homocysteine as a risk factor for atherosclerosis: A review. *Cardiovasc Pathol* 1997;6:1-9.
23. Lussier-Cacan S, Xhignesse M, Piolot A, Selhub J, Davignon J, Genest J Jr. Plasma total homocysteine in healthy subjects: sex-specific relation with biological traits. *Am J Clin Nutr* 1996; 64:587-593.
24. Clarke R, Woodhouse P, Ulvik A, Frost C, Sherliker P, Refsum H, Ueland PM, Khaw K- T. Variability and determinants of total homocysteine concentrations in plasma in an elderly population. *Clin Chem* 1998;44:102-107.
25. Jacques PF, Selhub J, Bostom AG, Wilson PWF, Rosenberg IH. The effect of folic acid fortification on plasma folate and total homocysteine concentrations. *N Engl J Med* 1 999;340: 1449-1454.
26. Lawrence JM, Pettitt DB, Watkins M, Umekubo MA. Trends in serum folate after food fortification. *Lancet* 1 999;354:91 5-916.
27. Fiskerstrand T, Refsum H, Kvalheim G, Ueland PM, Homocysteine and other thiols in plasma and urine: Automated determination and sample stability. *Clin Chem* 1 993;39:263-271.

Rev. (2011-11ax)cc