

Human High sensitivity IL-2 ELISA KIT

INTENDED USE

The Human high sensitivity IL-2 ELISA is to be used for the in-vitro quantitative determination interleukin-2 (IL-2) in human sera, plasmas, buffered solutions and cell culture media. The assay will recognize both natural and recombinant human IL-2. **This kit is intended for research use only and is not to be used in diagnostic procedures.**

PRINCIPLE OF THE METHOD

The IL-2 Kit is a solid phase sandwich Enzyme Linked-Immuno-Sorbent Assay (ELISA). A monoclonal antibody specific for IL-2 has been coated onto the wells of the microtiter strips provided. Samples, including standards of known IL-2 concentrations and unknowns are pipetted into these wells.

During the first incubation, the IL-2 antigen and a biotinylated monoclonal antibody specific for IL-2 are simultaneously incubated. After washing, the enzyme (streptavidin-peroxydase) is added. All the unbound enzyme is removed by washing and the first amplification step is performed by adding the Biotin-Tyramine reagent. Under the action of HRP, a biotin polymerisation reaction occurs in the region of the HRP linked to the detection antibody. After washing the second amplification step is performed and the polymerised biotin is revealed by a new streptavidin-HRP step. Finally after washing, the substrate is added. The intensity of this coloured product is directly proportional to the concentration of IL-2 present in the samples.

REAGENTS PROVIDED AND RECONSTITUTION

| REAGENTS (Store at 2-8°C) | COLOUR CODE | 1x96 wells Cat # 850.870.096 | 2x96 wells Cat # 850.870.192 | RECONSTITUTION |
|--------------------------------|-------------------|---------------------------------|---------------------------------|---|
| 96-wells microtiter plates | - | 1 | 2 | Ready-to-use. |
| Plastic cover | - | 2 | 4 | - |
| Standard : 60 pg/ml | Yellow | 2 vials | 4 vials | Reconstitute with the volume of standard diluent indicated on the vial. (See Reagents Preparation). |
| Standard Diluent Buffer | Black | 1 vial | 1 vial | (25 ml) 10X concentrate. Dilute in distilled Water. |
| Standard Diluent : human serum | Black | 1 vial | 2 vials | (7 ml) Ready-to-use. |
| Biotinylated anti-IL-2 | Red | 1 vial | 2 vials | (0.4 ml) Dilute in biotinylated antibody diluent. |
| Biotinylated Antibody Diluent | Red | 1 vial (7.5 ml) | 1 vial (13 ml) | Ready-to-use. |
| Streptavidin-HRP | | 2 vials | 4 vials | (5 µl) - Add 0.5 ml of HRP-Diluent before further dilutions. |
| Amplification Diluent | Brown & blue spot | 1 vial | 1 vial | (25 ml) Ready to use. |
| Amplifier* | Yellow | 1 vial | 2 vials | (200 µl) Dilute in Amplification buffer. |
| HRP Diluent | Red | 1 vial | 2 vials | (25 ml) Ready-to-use |
| Washing Buffer | White | 1 vial | 2 vials | (10 ml) 200X concentrate. Dilute in distilled Water . |
| Chromogen TMB : | Brown | 1 vial (11 ml) | 1 vial (24 ml) | Ready-to-use. |
| H2SO4 : Stop Reagent | Black | 1 vial | 2 vials | (11 ml) Ready-to-use. |

*Reagent contains ethyl alcohol.

MATERIAL REQUIRED BUT NOT PROVIDED

- Distilled water.
- Pipettes : 10 µl, 50 µl, 100 µl, 200 µl and 1000 µl.
- Vortex mixer and magnetic stirrer.
- Rotator set

SAFETY

- For research use only.
- The human blood components included in this kit have been tested and found non-reactive for HBsAg and anti-HIV. Nevertheless, no known method can offer complete assurance that human blood derivatives will not transmit hepatitis, AIDS or other infections. Therefore, handling of reagents, serum or plasma specimens should be in accordance with local safety procedures, e.g. CDC/NIH Health manual : "Biosafety in Microbiological and Biomedical Laboratories 4th Edition" 1999.
- Avoid any skin contact with H₂SO₄ and TMB. In case of contact, wash thoroughly

- Do not eat, drink, smoke or apply cosmetics where kit reagents are used.
- Do not pipette with mouth.

PROCEDURAL NOTES/LAB. QUALITY CONTROL

1. When not in use, kit components should be stored refrigerated as indicated on vials or bottles labels. All reagents should be warmed to room temperature before use. Lyophilized standards and controls should be discarded directly after resuspension and use.
2. Once the desired number of strips has been removed, immediately reseal the bag to protect the remaining strips from deterioration.
3. Cover or cap all reagents when not in use.
4. Do not mix or interchange reagents between different batches.
5. Do not use reagents beyond the expiration date of the kit.
6. Use a clean disposable plastic pipette tip for each reagent, standard, or specimen addition in order to avoid cross-contamination; for the dispensing of H₂SO₄ and substrate solution, avoid pipettes with metal parts.
7. Use a clean plastic container to prepare the washing solution.
8. Thoroughly mix the reagents and samples before use by agitation or swirling.
9. All residual washing liquid must be drained from the wells by efficient aspiration or by decantation followed by tapping the plate forcefully on absorbent paper. Never insert absorbent paper directly into the wells.
10. The TMB solution is light sensitive. Avoid prolonged exposure to light. Also, avoid contact of the TMB solution with metal to prevent colour development. Warning TMB is toxic avoid direct contact with hands. Dispose of properly.
11. When pipetting reagents, maintain a consistent order of addition from well-to-well. This will ensure equal incubation times for all wells.
12. Respect incubation times described in the assay procedure.

SPECIMEN COLLECTION, PROCESSING AND STORAGE

Cell culture supernatants- Remove particulates and aggregates by spinning at approximately 1000 x g for 10 min.

Serum–Avoid any unintentional stimulation of the cells by the procedure. Use pyrogen/endotoxin free collecting tubes. Serum should be removed rapidly and carefully from the red cells after clotting. For that, after clotting, centrifuge at approximately 1000 x g for 10 min and remove serum.

Plasma–EDTA, citrate and heparin plasma can be assayed. Spin samples at 1000 x g for 30 min to remove particulates. Harvest plasma.

Storage–If not analyzed shortly after collection, samples should be aliquoted (250-500µl) to avoid freeze-thaw cycles and stored frozen at –70°C. Avoid multiple freeze-thaw cycles of frozen specimens.

When possible, avoid use of badly hemolyzed or lipemic sera. If large amounts of particles are present, they should be removed prior to assay by centrifugation or filtration.

Recommendation : Do not thaw by heating at 37°C or 56°C. Thaw at room temperature and make sure that sample is completely thawed and homogeneous before assaying

PREPARATION OF REAGENTS

Standard buffer diluent 10X concentrate

Dilute 10 times with distilled water before use.

Standards

Depending on the type of samples you are assaying, the kit includes two standard diluents. Because biologicals fluids might contain proteases or cytokine-binding proteins that could modify the recognition of the cytokine you want to measure. You should reconstitute standard vials with the most appropriate standard diluent.

For serum and plasma samples use standard diluent human serum and for cells culture supernatants use Standard diluent buffer. Standard vials have to be reconstituted with the volume of appropriate standard diluent indicated on the vial. This reconstitution gives a stock solution of 60 pg/ml IL-2. Allow standard to stand for 5 minutes with gentle swirling prior to making dilutions. Serial dilutions of standard must be made before each assay and cannot be stored.

Samples


Normal sera and plasmas may be applied undiluted. Nevertheless, sera or plasmas from patients with various pathologies may be applied undiluted and diluted (to prevent too high concentrations). As IL-2 concentrations may vary considerably in cell supernatant samples, it is not easy to recommend a dilution factor. For example, unknown cell supernatant samples may also be tested undiluted and diluted.

Dilution of biotinylated anti-IL-2

Dilute the biotinylated anti-IL-2 with the biotinylated antibody diluent in a clean glass vial according to the number of wells to be used. See the next table for volumes to pipette. Extemporaneous preparations are recommended.

| Number of Wells | Biotinylated Antibody (µl) | Biotinylated Antibody Diluent (µl) |
|-----------------|----------------------------|------------------------------------|
| 16 | 40 | 1060 |
| 32 | 80 | 2120 |
| 48 | 120 | 3180 |
| 96 | 240 | 6360 |

Dilutions of Streptavidin-HRP – Step J and P :

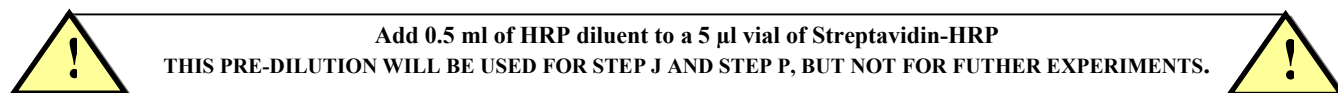

From the previous solution two others dilutions must be prepared: one for the step “J” and one for the step “P”.

Extemporaneous preparations are recommended. Following the number of wells to be used, further dilutions of Streptavidin-HRP should be made with HRP diluent in a clean glass vial: see hereafter the table for volumes to pipette.

| Number of Wells | Preparation of Streptavidin solution 1 – Step J | | Preparation of Streptavidin solution 2 – Step P | |
|-----------------|---|------------------------|---|------------------------|
| | Streptavidin-HRP (µl) | Strep-HRP Diluent (ml) | Streptavidin-HRP (µl) | Strep-HRP Diluent (ml) |
| 16 | 20 | 1.980 | 50 | 1.950 |
| 32 | 40 | 3.960 | 100 | 3.900 |
| 48 | 60 | 5.940 | 150 | 5.850 |
| 96 | 120 | 11.880 | 300 | 11.700 |

Dilution of Amplifier

Extemporaneous preparations are recommended. Following the number of wells to be used, further dilutions of amplifier should be made with amplification diluent in a clean glass vial: see hereafter the table for volumes to pipette.

| Number of Wells | Amplifier (µl) | Amplification Diluent (ml) |
|-----------------|----------------|----------------------------|
| 16 | 20 | 1.980 |
| 32 | 40 | 3.960 |
| 48 | 60 | 5.940 |
| 96 | 120 | 11.880 |

Washing Buffer 200X concentrate

Dilute 200 times in distilled water.

ASSAY METHOD

- a) Before use, mix all reagents thoroughly without making foam.
- b) Determine the number of microwell strips required to test the desired number of samples, plus appropriate number of wells needed for running blanks and standards. **Each sample, standard, blank should be assayed in duplicate.** Remove sufficient microwell strips from the pouch.
- c) Add 100 µl of appropriate standard diluent to standard wells B1, B2, C1, C2, D1, D2, E1, E2, F1, and F2. Reconstitute standard vial with the appropriate volume as described in the chapter reagents preparation. Pipet 200 µl of standard into wells A1 and A2 (see Plate Scheme below). Transfer 100 µl from A1 and A2 to B1 and B2 wells. Mix the contents by repeated aspirations and ejections. Take care not to scratch the inner surface of microwells. Repeat this procedure from the wells B1, B2 to wells C1, C2 and from wells C1, C2 to D1, D2 and so on creating two parallel rows of IL-2 standard dilutions ranging from 60 to 1.87 pg/ml. Discard 100 µl from the content of the last microwells used (F1, F2).
Alternatively these dilutions can be made in separate tubes and diluted standard pipetted directly into wells.
- d) Add 100 µl of appropriate standard diluent to the blank wells (G1-G2).
- e) Add 100 µl of sample to sample wells
- f) Preparation of biotinylated anti-IL-2: see reagents preparation.
- g) Add 50 µl of diluted biotinylated anti-IL-2 to all wells.
- h) Cover with a plate cover and incubate for 1 hour at room temperature (18°C - 22°C) with slow shaking
- i) Remove the cover and wash the plate as follows:
 - 1) aspirate the liquid from each well ;
 - 2) dispense 0.3 ml of washing solution into each well ;
 - 3) aspirate again the content of each well ;
 - 4) Repeat steps 2) and 3) two times.
- j) Prepare streptavidin-HRP solution 1 just before use: see reagents preparation.
- k) Distribute 100 µl of the streptavidin-HRP solution 1 to all wells, including blank wells.
- l) Cover and incubate 20 min at room temperature with slow shaking.
- m) Remove the cover and wash according to step i)
- n) Add 100 µl of amplifier dilution prepared according the instructions for the preparation of reagents and incubate 15 minutes with slow shaking.
- o) Remove the cover and wash according to step i).

All incubation steps except substrate TMB and acid are performed under slow shaking at room temperature

- p) Prepare streptavidin-HRP solution 2 just before use: see reagents preparation.
- q) Distribute 100 µl of the streptavidin-HRP solution 2 to all wells, including blank wells.
- r) Cover and incubate 20 min at room temperature with slow shaking.
- s) Remove the cover and wash according to step i).
- t) Pipette 100 µl of ready-to-use TMB substrate solution into all wells, including the blank wells and incubate in the dark for around 5 minutes at room temperature. Avoid direct exposure to light by wrapping the plate in aluminium foil.
- u) Incubation time of the substrate solution is usually determined by the ELISA reader performances: many ELISA readers record absorbance only up to 2.0 O.D. **The O.D. values of the plate should be monitored and the substrate reaction stopped before positive wells are no longer properly readable.**
- v) The enzyme-substrate reaction is stopped by quickly pipetting 100 µl of stop reagent (1N sulphuric acid) into each well, including the blank wells, to completely and uniformly inactivate the enzyme. Results must be read rapidly after the addition of stop reagent.
- w) Read absorbance of each well on a spectrophotometer using 450 nm as the primary wavelength and optionally 620 nm (610 nm to 650 nm is acceptable) as the reference wavelength.


Note: In case of incubation without shaking the O.D values may be lower than with shaking; in this case let the color develop longer in order to obtain correct OD values

SUGGESTED PLATE SCHEME

| | Standard Concentrations pg/mL | | Sample wells | | | | | | | | | |
|---|----------------------------------|-------|--------------|---|---|---|---|---|---|----|----|----|
| | 1 | 2 | 3 | 4 | 5 | 6 | 7 | 8 | 9 | 10 | 11 | 12 |
| A | 60 | 60 | | | | | | | | | | |
| B | 30 | 30 | | | | | | | | | | |
| C | 15 | 15 | | | | | | | | | | |
| D | 7.5 | 7.5 | | | | | | | | | | |
| E | 3.75 | 3.75 | | | | | | | | | | |
| F | 1.87 | 1.87 | | | | | | | | | | |
| G | Blank | Blank | | | | | | | | | | |
| H | | | | | | | | | | | | |

DATA ANALYSIS

Generate a linear standard curve by plotting the average absorbance on the vertical axis versus the corresponding IL-2 standard concentration on the horizontal axis. The amount of IL-2 in each sample is determined by extrapolating OD values to IL-2 concentrations using the standard curve.

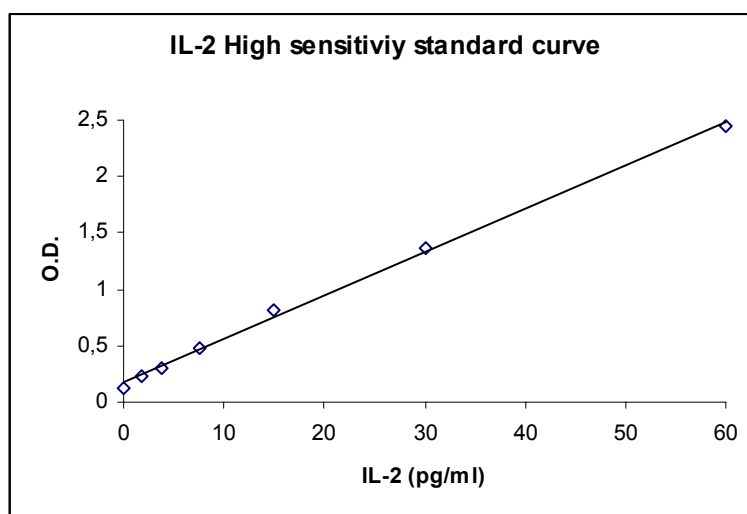

Typical IL-2 High sensitivity standard curve ranging from 1.87 to 60 pg/mL

LIMITATIONS OF THE PROCEDURE

Do not extrapolate the standard curve beyond the 60 pg/ml standard curve point. The dose-response is non-linear in this region and therefore it is difficult to obtain accurate concentration. Concentrated samples (> 60 pg/ml) have to be diluted with standard diluent or with your own sample buffer. During analysis, multiply results by the appropriate dilution factor.
The influence of various drugs, aberrant sera (hemolyzed, hyperlipidemic, jaundiced...) has not been investigated. The rate of degradation of native IL-2 in various matrices has not been investigated.

PERFORMANCES AND CHARACTERISTICS

Sensitivity

This has been determined by adding 3 standard deviations to the mean concentration of 40 zeros. The minimum detectable dose of IL-2 was less than 0.97 pg/ml.

Specificity

Ten specificities were tested with concentrations higher than IL-2 curve concentrations. No cross reaction was observed for concentrations ranging from 300 pg/ml to 18.75 pg/ml for IL-1 α and β , IL-8, IL-10, TNF- α , GrB, FasL, TRAIL, CD116 and Gp-130.

Spike - Recovery

The spike recovery was evaluated by spiking three levels of IL-2 into three different pooled human sera and one cell culture medium. Recovery was evaluated with one test. The recovery in pooled human sera ranged from 89 to 111% with an average of 100%. In cell culture medium, the recovery is 95%.

Precision

Four sera and one cell culture medium samples spiked with various concentrations of IL-2 were tested for repeatability and reproductibility. Each assay was carried out with 3 duplicates of each sample. Five independent assays were performed. The intra-assay and inter-assay coefficient of variation were calculated to be 4.2% and 9.0% respectively. Note that precision in plasma sample is lower than precision in sera or culture medium (data not shown).

Linearity of dilution

Two pooled human sera, one plasma and one cell culture medium containing different concentrations of IL-2 were serially diluted in standard buffer diluent. Linearity was evaluated on 4 dilutions. The linear regression of samples versus the expected concentration yielded a quote slope of 0.996.

Expected values

16 sera from healthy individual donors were tested undiluted in duplicates. 13 sera were negatives and 3 are slightly positives: 1.33, 1.59 and 2.36 pg/ml.

>ε-----

Study : _____

Specificity : _____

Date : _____

Plate n°: _____

| | 1 | 2 | 3 | 4 | 5 | 6 | 7 | 8 | 9 | 10 | 11 | 12 |
|---|---|---|---|---|---|---|---|---|---|----|----|----|
| A | | | | | | | | | | | | |
| B | | | | | | | | | | | | |
| C | | | | | | | | | | | | |
| D | | | | | | | | | | | | |
| E | | | | | | | | | | | | |
| F | | | | | | | | | | | | |
| G | | | | | | | | | | | | |
| H | | | | | | | | | | | | |

ASSAY PROCEDURE SUMMARY :

Human High Sensitivity IL-2 ELISA Kit

Total procedure length: 2H00

| | |
|--|---|
| Specific molecule detection steps | <p>Add 100µl of sample or diluted standard</p> <p style="text-align: center;">⇓</p> <p>Add 50µl of diluted biotinylated detection antibody to all wells</p> <p style="text-align: center;">⇓</p> <p>Incubate 1 hour at room temperature with slow shaking</p> <p style="text-align: center;">⇓</p> |
| Streptavidin –HRP / Amplification steps | <p style="text-align: center;">⇓</p> <p>Wash three times</p> <p style="text-align: center;">⇓</p> <p>Add 100µl of streptavidin-HRP (solution 1) to all wells</p> <p style="text-align: center;">⇓</p> <p>Incubate 20 min at room temperature with slow shaking</p> <p style="text-align: center;">⇓</p> <p>Wash three times</p> <p style="text-align: center;">⇓</p> <p>Add 100 µl amplifier to all wells</p> <p style="text-align: center;">⇓</p> <p>Incubate 15 min at room temperature with slow shaking</p> <p style="text-align: center;">⇓</p> <p>Wash three times</p> <p style="text-align: center;">⇓</p> <p>Add 100µl of streptavidin-HRP (solution 2) to all wells</p> <p style="text-align: center;">⇓</p> <p>Incubate 20 min at room temperature with slow shaking</p> <p style="text-align: center;">⇓</p> |
| Revelation and reading steps | <p style="text-align: center;">⇓</p> <p>Wash three times</p> <p style="text-align: center;">⇓</p> <p>Add 100 µl of ready-to-use TMB Protect from light. Let the color develop for around 5 min.</p> <p style="text-align: center;">⇓</p> <p>Add 100 H2SO4</p> <p style="text-align: center;">⇓</p> <p>Read Absorbance at 450 nm</p> |