

Leptin Human E. coli Tag free

Product Data Sheet

Type: Recombinant

Tag: Tagless

Source: E. coli

Species: Human

Other names: Obesity factor, Obese protein, LEP, OB, OBS

Cat. No.:

RD172001100

(0.1 mg)

Description

Total 147 AA. MW: 16 kDa (calculated). UniProtKB acc.no. P41159.

Introduction to the Molecule

Leptin, the product of the *ob* (obese) gene, is a single-chain 16 kDa proteohormone consisting of 146 amino acid residues.

Leptin is produced by differentiated adipocytes, although production have been demonstrated in other tissues, such as fundus of the stomach, the skeletal muscle, the liver, and the placenta. Leptin is considered to play an important role in appetite control, fat metabolism and body weight regulation. It targets the central nervous system, in particular the hypothalamus, suppressing food intake and stimulating energy expenditure. In humans, leptin levels correlate with body mass index (BMI) and percentage body fat, and are elevated even in obese individuals. Leptin has a dual action; it decreases the appetite and increases energy consumption, causing more fat to be burned.

Research topic

Diabetology - Other Relevant Products, Energy metabolism and body weight regulation, Reproduction

Amino Acid Sequence

AVPIQKVQDD TKTLIKTIVT RINDISHTQS VSSKQKVTGL DFIPGLHPIL TLSKMDQTLA VYQQILTSMR SRNVIQISND
LENLRDLLHV LAFSKSCHLP WASGLETLDS LGGVLEASGY STEVVALSRL QGSLQDMLWQ LDLSPGC


Source

E. coli

Purity

>95%

SDS-PAGE gel

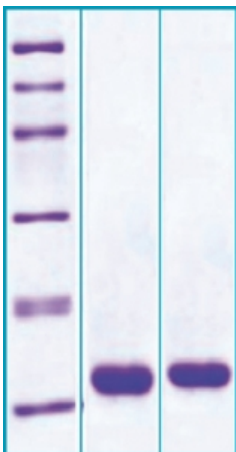

12% SDS-PAGE separation of Human Leptin

1. M.W. marker - 14, 21, 31, 45, 66, 97 kDa

2. reduced and heated sample, 7 µg/lane

3. non-reduced and non-heated sample, 7 µg/lane

Biological Activity

Biological Activity is evidenced by inducing proliferation of BAF/3 cells stably transfected with the long form of human leptin receptor.

Endotoxin

< 1.0 EU/ug

Formulation

Filtered (0.4 µm) and lyophilized from 0.5 mg/ml in 5mM NaHCO₃, pH= 7,4

Reconstitution

Add deionized water to prepare a working stock solution of approximately 0.5 mg/mL and let the lyophilized pellet dissolve completely. Filter sterilize your culture media/working solutions containing this non-sterile product before using in cell culture.

Shipping

At ambient temperature. Upon receipt, store the product at the temperature recommended below.

Storage, Stability/Shelf Life

Store lyophilized protein at -80°C. Lyophilized protein remains stable until the expiry date when stored at -80°C. Aliquot reconstituted protein to avoid repeated freezing/thawing cycles and store at -80°C for long term storage. Reconstituted protein can be stored at 4°C for a week.

Quality Control Test

BCA to determine quantity of the protein.

SDS PAGE to determine purity of the protein.

LAL to determine quantity of endotoxin.

Applications

ELISA, Western blotting

Note

This product is intended for research use only.

References to this Product

- Stejskal D, Karpisek M . *Adipocyte fatty acid binding protein in a Caucasian population: a new marker of metabolic syndrome?*. [Eur J Clin Invest](#) . Sep;36(9):621-5 (2006)
- Nonaka S, Hashizume T, Yamashita T . *Effects of leptin and leptin peptide amide on the release of luteinizing hormone, growth hormone and prolactin from cultured porcine anterior pituitary cells*. [Animal Science Journal](#) . 77, 44-52 (2006)
- Johnstone AM, Murison SD, Duncan JS, Rance KA, Speakman, JR . *Factors influencing variation in basal metabolic rate include fat-free mass, fat mass, age, and circulating thyroxine but not sex, circulating leptin, or triiodothyronine*. [Am J Clin Nutr](#) . Nov;82(5):941-8 (2005)
- Matarese G, Carrieri PB, La Cava A, Perna F, Sanna V, De Rosa V, Aufiero D, Fontana S, Zappacosta S . *Leptin increase in multiple sclerosis associates with reduced number of CD4(+)/CD25+ regulatory T cells*. [Proc Natl Acad Sci U S A](#) . Apr 5;102(14):5150-5 (2005)
- Lee MJ, Fried SK . *Multilevel regulation of leptin storage, turnover, and secretion by feeding and insulin in rat adipose tissue*. [J Lipid Res](#) . Sep;47(9):1984-93 (2006)
- Kim EJ, Kim YK, Kim MK, Kim S, Kim JY, Lee DH, Chung JH. *UV-induced inhibition of adipokine production in subcutaneous fat aggravates dermal matrix degradation in human skin*. [Sci Rep](#). 2016;6:25616

HEADQUARTERS: BioVendor Laboratorní medicína, a.s.	Karasek 1767/1	621 00 Brno CZECH REPUBLIC	Phone: +420-549-124-185 Fax: +420-549-211-460	E-mail: info@biovendor.com sales@biovendor.com Web: www.biovendor.com
AUSTRIA: BioVendor GesmbH	Nußdorfer Straße 20/10	1090 Vienna AUSTRIA	Phone: +43-1-89090-25 Fax: +43-1-89090-2515	E-mail: infoAustria@biovendor.com
GERMANY, SWITZERLAND: BioVendor GmbH	Otto-Hahn-Straße 16	34123 Kassel GERMANY	Phone: +49-6221-433-9100 Fax: +49-6221-433-9111	E-mail: infoEU@biovendor.com
USA, CANADA AND MEXICO: BioVendor LLC	128 Bingham Rd. Suite 1300	Asheville, NC 28806 USA	Phone: +1-828-575-9250 +1-800-404-7807 Fax: +1-828-575-9251	E-mail: infoUSA@biovendor.com