

## Connective Tissue Growth Factor Human E. coli

### Product Data Sheet

<b>Type:</b> Recombinant	<b>Cat. No.:</b>	
<b>Source:</b> E. coli	RD172035025	(0.025 mg)
<b>Species:</b> Human	RD172035100	(0.1 mg)
<b>Other names:</b> Hypertrophic chondrocyte-specific protein 24, Insulin-like growth factor-binding protein 8, IGF-binding protein 8, IGFBP-8, IBP-8, CCN family member 2, CTGF, CCN2, HCS24		

### Description

Total 346 AA. MW: 38.3 kDa (calculated). UniProtKB acc.no. P29279. N-Terminal His-tag and Xa - cleavage site, 23 extra AA (highlighted).

### Introduction to the Molecule

Connective Tissue Growth Factor belongs to the CCN family of proteins. The CCN family presently consists of six members in human also known as: Cyr61 (Cystein rich 61), CTGF (Connective Tissue Growth Factor), Nov (Nephroblastoma Overexpressed gene), WISP-1, 2 and 3 (Wnt-1 Induced Secreted Proteins). The CCN genes encode secreted proteins associated with the Extracellular Matrix (ECM) and cell membrane.

CCN proteins are matricellular proteins which are involved in the regulation of various cellular functions including: proliferation, differentiation, survival, adhesion and migration. They are expressed in derivatives of the three embryonic sheets and are implicated in the development of kidney, nervous system, muscle, bone marrow, cartilage and bone. During adulthood, they are implicated in wound healing, bone fracture repair, and pathologies such as: fibrosis, vascular ailments and tumourigenesis. Full length secreted CCN proteins can show an antiproliferative activity, whereas truncated isoforms are likely to stimulate proliferation and behave as oncogenes.

The full-length protein consists of four modules:

- Module I shares partial identity with the N-terminal part of the Insulin-like Growth Factor Binding Proteins (IGFBPs).
- Module II includes a stretch of 70 amino acid residues - which shares sequence identity with the Von Willebrand Factor Type C repeat (VWC).
- Module III contains sequences sharing identity with the Thrombospondin type 1 repeat (TSP1) (WSXCSXXCG), which is thought to be implicated in the binding of sulfated glycoconjugates and to be important for cell adhesion.
- Module IV, also designated CT, is encoded by exon 5. It is the least conserved one of the four domains at the level of nucleotide sequence, but it appears to be critical for several of the biological functions attributed to the CCN proteins. Module IV resembles the CT domain of several extracellular protein including, Von Willebrand's factor and mucins. Sequence similarities to heparin-binding motifs are also found within this domain.

Proteolysis of the secreted full-length CCN proteins that has been reported in the case of CCN2 and CCN3 might result in the production of CCN-derived peptides with high affinity for ligands that full-length CCN proteins bind only poorly. Amino-truncated CCN2 isoforms were biologically active whereas no specific biological activity has been attributed to the truncated CCN3. Although the molecular processes underlying the production of these secreted isoforms is presently unknown, it is important to note that proteolysis occur at the same amino acid residues in both CCN2 and CCN3. An elevated expression of CCN2 has also been detected by Northern blotting in human invasive mammary ductal carcinomas, dermatofibromas, pyogenic granuloma, endothelial cells of angioliomas and angioleiomyomas, and in pancreatic tumours. A study performed with chondrosarcomas representative of various histological grades established that CCN2 expression was closely correlated with increasing levels of malignancy. In agreement with CCN2 playing a role in brain tumour angiogenesis, immunocytochemistry studies indicated that both glioblastoma tumour cells and proliferating endothelial cells stained positive for CCN2. In astrocytomas, CCN2 expression was particularly elevated in high-grade tumours, with a marked effect of CCN2 on cell proliferation. Downregulation of CCN2 expression in these cells was associated with a growth arrest at the G1/S transition while over-expression of CCN2 induced a two-fold increase of the number of cells in the G1 phase. Gene profiling analysis allowed to identify a set of about 50 genes whose expression might account for the proliferative activity of CCN2 in these cells. CCN2 was seen in a higher proportion of mononuclear cells of patients with acute lymphoblastic leukemia.

## Research topic

Others

## Amino Acid Sequence

MGHHHHHHHH HHSSGHIEGR HMRQNCSGPC RCPDEFAPRC PAGVSLVLDG CGCCRVCAKQ LGELCTERDP CDPHKGLFCD  
FGSPANRKIG VCTAKDGAPC IFGGTVYRSG ESFQSSCKYQ CTCLDGAVGC MPLCSMDVRL PSPDCPFPRR VKLPGKCCEE  
WVCDEPKDQT VVGPALAAYR LEDTFGPDPT MIRANCLVQT TEWSACSKTC GMGISTRVTN DNASCRLEKQ SRLCMVRPCE  
ADLEENIKKG KKCIRTPKIS KPIKFELSGC TSMKTYRAKF CGVCTDGRCC TPHRTTTLPV EFKCPDGEVM KKNMMFIKTC  
ACHYNCPGDN DIFESLYYRK MYGDMA


## Source

E. coli

## Purity

Purity as determined by densitometric image analysis: > 95 %

## SDS-PAGE gel

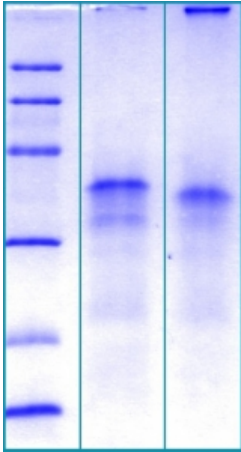

12% SDS-PAGE separation of Human CTGF

1. M.W. marker - 14, 21, 31, 45, 66, 97 kDa

2. reduced and heated sample, 2.5µg/lane

3. non-reduced and non-heated sample, 2.5µg/lane

## Endotoxin

< 1.0 EU/ug

## Formulation

Filtered (0.4 µm) and lyophilized in 0.5 mg/mL in 0.05 M Acetate buffer pH=4.0; 4%mannitol, 1% sucrose

## Reconstitution

Add 0.1M Acetate buffer pH=4.0 to prepare a working stock solution of 0.5 mg/mL and let the lyophilized pellet dissolve completely at 37°C. For conversion into higher pH value, we recommend intensive dilution by relevant buffer to a concentration of 10µg/mL. In higher concentrations the solubility of this antigen is limited. Filter sterilize your culture media/working solutions containing this non-sterile product before using in cell culture.

## Shipping

On ice. Upon receipt, store the product at the temperature recommended below.

## Storage, Stability/Shelf Life

Store lyophilized protein at -80°C. Lyophilized protein remains stable until the expiry date when stored at -80°C. Aliquot reconstituted protein to avoid repeated freezing/thawing cycles and store at -80°C for long term storage. Reconstituted protein can be stored at 4°C for a limited period of time; it does not show any change after one week at 4°C.

## Quality Control Test

BCA to determine quantity of the protein.

SDS PAGE to determine purity of the protein.

LAL to determine quantity of endotoxin.

## Applications

Cell culture and/or animal studies, Western blotting

## Note

This product is intended for research use only. The Certificate of Analysis is available on [www.biovendor.com](http://www.biovendor.com)

## References to this Product

- Yoshida K, Munakata H . *Connective tissue growth factor binds to fibronectin through the type I repeat modules and enhances the affinity of fibronectin to fibrin.* [Biochim Biophys Acta](#) . Apr;1770(4):672-80 (2007)
- Gressner AM, Yagmur E, Lahme B, Gressner O, Stanzel S . *Connective tissue growth factor in serum as a new candidate test for assessment of hepatic fibrosis.* [Clin Chem](#) . Sep;52(9):1815-7 (2006)
- Hayata N, Fujio Y, Yamamoto Y, Iwakura T, Obana M, Takai M, Mohri T, Nonen S, Maeda M, Azuma J. *Connective tissue growth factor induces cardiac hypertrophy through Akt signaling.* *Biochem Biophys Res Commun.* 2008 Mar 27;
- Lee CH, Shah B, Moiola EK, Mao JJ. *CTGF directs fibroblast differentiation from human mesenchymal stem/stromal cells and defines connective tissue healing in a rodent injury model.* *J Clin Invest.* 2010 Sep;120 (9):3340-9
- Hu R, Ling W, Xu W, Han D. *Fibroblast-like cells differentiated from adipose-derived mesenchymal stem cells for vocal fold wound healing.* *PLoS One.* 2014;9 (3):e92676
- Lee MS, Ghim J, Kim SJ, Yun YS, Yoo SA, Suh PG, Kim WU, Ryu SH. *Functional interaction between CTGF and FPRL1 regulates VEGF-A-induced angiogenesis.* *Cell Signal.* 2015 Jul;27 (7):1439-48
- Aoyama E, Hattori T, Hoshijima M, Araki D, Nishida T, Kubota S, Takigawa M. *N-terminal domains of CCN family 2/connective tissue growth factor bind to aggrecan.* *Biochem J.* 2009 Jun 15;420 (3):413-20
- Aoyama E, Hattori T, Hoshijima M, Araki D, Nishida T, Kubota S, Takigawa M. *N-terminal domains of CCN family 2/connective tissue growth factor bind to aggrecan.* *Biochem J.* 2009 Jun 15;420 (3):413-20

<b>HEADQUARTERS:</b> BioVendor Laboratorní medicina, a.s.	Karasek 1767/1	621 00 Brno CZECH REPUBLIC	Phone: +420-549-124-185 Fax: +420-549-211-460	E-mail: info@biovendor.com sales@biovendor.com Web: www.biovendor.com
AUSTRIA: BioVendor GesmbH	Nußdorfer Straße 20/10	1090 Vienna AUSTRIA	Phone: +43-1-89090-25 Fax: +43-1-89090-2515	E-mail: infoAustria@biovendor.com
GERMANY, SWITZERLAND: BioVendor GmbH	Otto-Hahn-Straße 16	34123 Kassel GERMANY	Phone: +49-6221-433-9100 Fax: +49-6221-433-9111	E-mail: infoEU@biovendor.com
USA, CANADA AND MEXICO: BioVendor LLC	128 Bingham Rd. Suite 1300	Asheville, NC 28806 USA	Phone: +1-828-575-9250 +1-800-404-7807 Fax: +1-828-575-9251	E-mail: infoUSA@biovendor.com