

Resistin Mouse E. coli

Product Data Sheet

Type: Recombinant

Source: E. coli

Species: Mouse

Cat. No.:

RD272016300

(0.1 mg)

Other names: Cysteine-rich secreted protein FIZZ3, Adipose tissue-specific secretory factor, ADSF, C/EBP-epsilon-regulated myeloid-specific secreted cysteine-rich protein, Cysteine-rich secreted protein A12-alpha-like 2, RETN, FIZZ3, HXCP1, RSTN, UNQ407/PRO1199

Description

Total 103 AA. MW: 11.3 kDa (calculated). UniProtKB acc. No. Q99P87 (Ser21-Ser114). N-terminal His-tag (9 extra AA). Protein identity confirmed by LC-MS/MS.

Introduction to the Molecule

Resistin, a product of the RSTN gene, is a peptide hormone belonging to the class of cysteine-rich secreted proteins which is termed the RELM family, and is also described as ADSF (Adipose Tissue-Specific Secretory Factor) and FIZZ3 (Found in Inflammatory Zone). Human resistin contains 108 amino acids as a prepeptide, and its hydrophobic signal peptide is cleaved before its secretion. Resistin circulates in human blood as a dimeric protein consisting of two 92 amino acid polypeptides, which are disulfide-linked via Cys26. Resistin may be an important link between obesity and insulin resistance. Mouse resistin, specifically produced and secreted by adipocyte, acts on skeletal muscle myocytes, hepatocytes and adipocytes themselves so that it reduces their sensitivity to insulin. Stepan et al. have suggested that resistin suppresses the ability of insulin to stimulate glucose uptake. They have also suggested that resistin is present at elevated levels in blood of obese mice, and is down regulated by fasting and antidiabetic drugs. Way et al., on the other hand, have found that resistin expression is severely suppressed in obesity and is stimulated by several antidiabetic drugs. Other studies have shown that mouse resistin increases during the differentiation of adipocytes, but it also seems to inhibit adipogenesis. In contrast, the human adipogenic differentiation is likely to be associated with a down regulation of resistin gene expression. Recent studies have shown that human resistin is expressed also in macrophages and may be a novel link between inflammation and insulin resistance.

Research topic

Animal studies, Diabetology - Other Relevant Products, Energy metabolism and body weight regulation

Amino Acid Sequence

MKHHHHHHAS SMP LCPIDEA IDKKIKQDFN SLFPNAIKNI GLNCWTVSSR GKLASCPEGT AVLSCSCGSA CGSWDIREEK
VCHCQCARID WTAARCKLQ VAS


Source

E. coli

Purity

Purity as determined by densitometric image analysis: >95%

SDS-PAGE gel

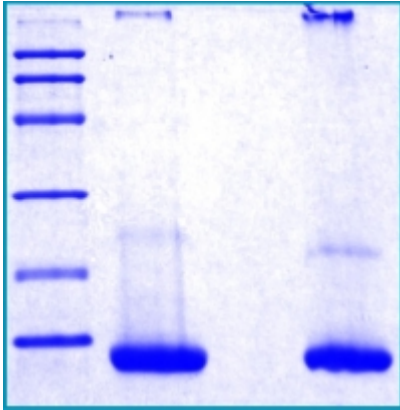

14% SDS-PAGE separation of Mouse Resistin
1. M.W. marker - 14, 21, 31, 45, 66, 97 kDa
2. reduced and heated sample, 2.5µg/lane
3. non-reduced and non-heated sample, 2.5µg/lane

Endotoxin

<0.1 EU/µg

Formulation

Filtered (0.4 µm) and lyophilized in 0.5 mg/mL in 0.05 M Acetate buffer pH=4.0

Reconstitution

Add 0.1M Acetate buffer pH4 to prepare a working stock solution of approximately 0.5 mg/mL and let the lyophilized pellet dissolve completely. For conversion into higher pH value, we recommend intensive dilution by relevant buffer to a concentration of 10µg/mL. In higher concentrations the solubility of this antigen is limited. Filter sterilize your culture media/working solutions containing this non-sterile product before using in cell culture.

Shipping

At ambient temperature. Upon receipt, store the product at the temperature recommended below.

Storage, Stability/Shelf Life

Store lyophilized protein at -80°C. Lyophilized protein remains stable until the expiry date when stored at -80°C. Aliquot reconstituted protein to avoid repeated freezing/thawing cycles and store at -80°C for long term storage. Reconstituted protein can be stored at 4°C for a limited period of time; it does not show any change after one week at 4°C.

Quality Control Test

BCA to determine quantity of the protein.
SDS PAGE to determine purity of the protein.
LAL TEST to determine quantity of endotoxin.

Applications

ELISA, Western blotting

Note

This product is intended for research use only. The recombinant protein is for research use only.

References to this Product

- Pravenec M, Kazdova L, Cahova M, Landa V, Zidek V, Mlejnek P, Simakova M, Wang J, Qi N, Kurtz TW . *Fat-specific transgenic expression of resistin in the spontaneously hypertensive rat impairs fatty acid re-esterification.* [Int J Obes \(Lond\)](#) . Jul;30(7):1157-9 (2006)

HEADQUARTERS:				
BioVendor Laboratorní medicína, a.s.	Karasek 1767/1	621 00 Brno CZECH REPUBLIC	Phone: +420-549-124-185 Fax: +420-549-211-460	E-mail: info@biovendor.com sales@biovendor.com Web: www.biovendor.com
AUSTRIA: BioVendor GesmbH	Nußdorfer Straße 20/10	1090 Vienna AUSTRIA	Phone: +43-1-89090-25 Fax: +43-1-89090-2515	E-mail: infoAustria@biovendor.com
GERMANY, SWITZERLAND: BioVendor GmbH	Otto-Hahn-Straße 16	34123 Kassel GERMANY	Phone: +49-6221-433-9100 Fax: +49-6221-433-9111	E-mail: infoEU@biovendor.com
USA, CANADA AND MEXICO: BioVendor LLC	128 Bingham Rd. Suite 1300	Asheville, NC 28806 USA	Phone: +1-828-575-9250 +1-800-404-7807 Fax: +1-828-575-9251	E-mail: infoUSA@biovendor.com